

BioTemps® *NOW*

"The Temporary You Admire"



Affordable PMMA material designed for **esthetic, in-office milled temporaries**


- Indicated for temporary full-coverage anterior and posterior crowns
- Mill in under 10 minutes*
- Rapid patient delivery with no oven wait
- Low water sorption
- Exclusively available for Glidewell milling systems



QUICK
MILL TIME



RETAIN GINGIVAL
POSITION

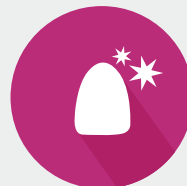

PROVISIONAL ANTERIOR
AND POSTERIOR CROWNS



PATIENT-PLEASING
ESTHETICS

For more information
glidewelldental.com ■ 800-854-7256



**GLIDEWELL
DENTAL**

BENEFITS

- Use BioTemps® NOW to let patients preview the final results
- Up to six months of wear time
- Shade stability long after delivery
- Opalescence and fluorescence that mimics natural teeth
- Excellent abrasion resistance
- Cost effective in-office milling of temporaries



Sold in 5-block packs

READY TO DELIVER IN UNDER 20 MINUTES



MILL

+



FINISH/POLISH

+



FURNACE
(No oven required)

+



CEMENT

=



TOTAL TIME

SPECIFICATIONS

MATERIAL

- PMMA
- Flexural Strength ≥ 50 MPa

INDICATIONS

- Temporary full-coverage anterior and posterior crowns

IDEAL FOR


- Provisional restorations

SHADES


- Available in six shades that correspond to the VITA® classical shade system and one bleach shade: A1, A2, A3, A3.5, B1, C2, BL

EASY PREPARATION

POSTERIOR

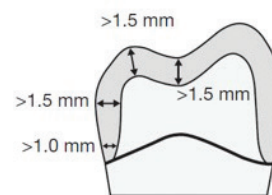

ANTERIOR

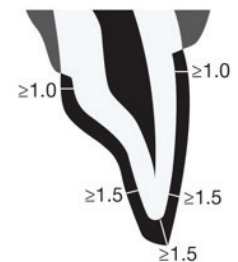

- Wall thickness should be at least 1.5 mm, 1.0 mm at margins.
- Prepare margins with a deep chamfer or a rounded shoulder.

*Milling times vary based on prescription and mill used. Average milling time is shown.

VITA is a registered trademark of VITA Zahnfabrik.

MKT 10465_3.0 PK-281915-092018

For more information
glidewelldental.com ■ 800-854-7256

